



Trolley for patient transport

# WP-09



## Trolley for patient transport

# WP-09

### Purpose

The trolley WP-09 is used for transporting patients between clinical areas. There is a possibility of monitoring a patient with a C-arm.



### Standard construction

- The trolley pallet is X-ray translucent. Thanks to a tray for X-ray cassette (option) which can be fixed under the trolley pallet it is possible to take X-ray pictures on the whole pallet length
- The trolley is made of steel sections covered with epoxy polyester powder paint in RAL 9002 colour resistant to scratches and disinfectant agents
- Back rest segment in the trolley and the thigh rest segment in the version of the trolley with 4 section pallet are adjusted with the use of gas springs with blockade
- The shank rest segment adjustment in the version of the trolley with 4 section pallet is carried out with the use of ratchet mechanism
- Base of the trolley is protected with ABS covers in RAL 9002 colour. It is possible to place an oxygen bottle on the trolley base
- The trolley pallet consists of two segments with mattress and basket designed for storing the patient belongings. Bumpers fixed in corners of the pallet protect it against a damage

### Trolley features

- Height adjustment of the pallet and longitudinal tilts are carried out with the use of hydraulic actuators adjusted by foot levers
- The height of the trolley pallet can be adjusted to the height of a hospital bed or an operating table
- The WP-09 trolley is equipped with four wheels with a central blockade. There is a possibility of fixing the fifth wheel. One of the wheels is a directional wheel and enables moving the trolley forward
- Big diameter of wheels (200 mm) facilitates manoeuvring the trolley
- Head board is removable

### Additional equipment

- Side rails
- Drip bottle holders
- Belt for fixing the patient
- Tray for X-ray cassette
- Fifth wheel



Movable tray placed under the trolley pallet for X-ray cassette WW-24.0 makes it possible to take X-ray pictures on the whole pallet length



Belt WW-12.0 for fixing the patient



Varnished side rails PB-15.5



Fifth wheel WW-25.0



Chromium plated side rails PB-15.6



Drip bottle holder WK-06.1



Drip bottle holder WK-16.0

## Technical data

Total length	2100 mm
Total width (without side rails)	800 mm
Total width (with side rails)	860 mm
Mattress dimensions	1920 x 675 x 60 mm
Pallet height adjustment	560 to 890 mm
Back rest inclination	0° to 70°
Thigh rest inclination (Version of the trolley with 4 section pallet)	30°
Shank rest inclination (Version of the trolley with 4 section pallet)	22°
Trendelenburg position	15°
Reverse Trendelenburg position	15°
Wheels diameter	200 mm
Permissible load (static)	250 kg
Permissible load (dynamic)	220 kg

**Famed Żywiec Sp. z o.o.**  
ul. Fabryczna 1  
34-300 Żywiec, Poland

**Head Office:**  
tel.: +48 33 866 62 00  
fax: +48 33 475 58 90



**Export Department:**  
tel.: +48 33 866 62 58  
tel.: +48 33 866 62 65  
tel.: +48 33 866 62 66  
fax: +48 33 861 46 78

expo.famed@famed.com.pl  
www.famed.com.pl




Functions of the trolley adjustable by a foot-operated pump

560 to 890 mm

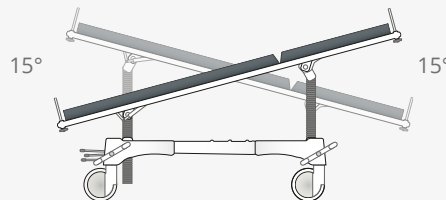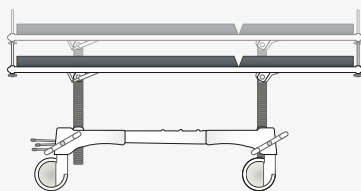

Functions of the trolley adjustable manually

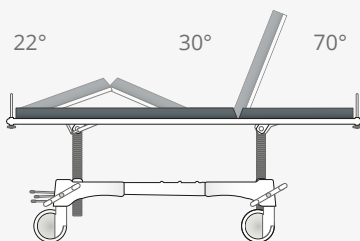

The product meets the requirements of the European Directive MDD 93/42/EEC concerning safety of medical devices

Edition 05/2015/01. Famed Żywiec reserves the right to make modifications in the products in the course of technical progress. All illustrations and photos, used in this material, are for illustration and may not accurately depict the actual product.

